



Non-Contact Tonometer

Kowa

KT-800




Technology for Life Science





Technology for Life Science

Simple operation with 3D auto alignment

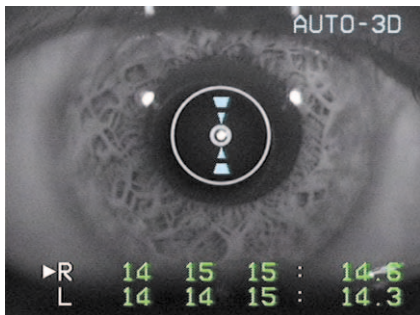

Special Features

3D auto alignment function

Up and down, right and left, and forward and backward positioning can be automatically performed only by placing the luminescent spot inside the circle.

Soft airflow

Improved softness of the air flow. Measurement can be done with the soft air feeling even from the first measurement.

No measurement range switching necessary

Smooth and easy measurement is possible without the range switching operation (30/60mmHg) as before.

Insufficient eye opening sensor function

Measurement is not processed if the eye is not opened sufficiently.

Easy paper replacement

Pursued usability and speedy printing.

Electrically operated chin rest

Height of the chin rest can be adjusted with the buttons at hand.



Specifications

| | |
|-----------------------------|---|
| Range of measurement | 0 - 60mmHg. (no range switching) |
| Unit of measurement | 1mmHg |
| Distance of measurement | 11mm |
| Monitor | 5.6 inch color TFT |
| Printer | 58mm width, thermal line printer |
| Frame moving range | Back and forth 40mm, left/right 86mm |
| Head moving range | Up/down 28mm |
| Auto alignment moving range | Up/down \pm 3mm, left/right \pm 3mm, back and forth \pm 2mm |
| Power supply | Input : AC100V - AC230V 50/60 Hz Power consumption : 60VA (normal)/ 110 VA (maximum) |
| Dimensions | 274 (W) x 457 (D) x 458 (H) mm |
| Weight | 18kg |

Distribution name : Kowa KT-800



Kowa Kowa Company, Ltd.

World Sales Headquarters

4-14, Nihonbashi-honcho 3-chome, Chuo-ku, Tokyo 103-8433 Japan
Phone: 81-3-3279-7331 Facsimile: 81-3-5255-7516
URL: <http://www.kowa.co.jp/e-life/>


Hamamatsu Factory

3-1, Shinmiyakoda 1-chome, Hamamatsu City,
Shizuoka Pref., 431-2103 Japan

Kowa Europe GmbH *Kowa Optimed, Inc.*

Immermannstrasse 43B
40210 Düsseldorf F.R.Germany
Phone : 49-211-1793540
Facsimile : 49-211-161952
URL : <http://kowa-europe.com/>

20001 S. Vermont Ave. Torrance,
CA 90502, U.S.A.
Phone : 1-310-327-1913
Facsimile : 1-310-327-4177
URL : <http://medical.kowa.com/>



This catalog was made with 100% recycled paper Printed in Japan 05971000MA

